




TECHNOCRATIC SERVICES

ISO 9001:2015 Certified Company

Laboratory Microbiology Thermostatic Heating Incubator

Description	Image
<p><u>Laboratory Microbiology Thermostatic Heating Incubator</u></p> <p><u>Feature:</u></p> <p>1.Precision LCD microcomputer intelligent PID control system, with temperature control ,timing and over-temperature alarm ,with RS485 computer interface.</p> <p>2. With standard international brands cooling fan, maintenance-free, durability. Perfect combination of duct design to ensure temperature uniformity for the working chamber.</p> <p>3.Mirror stainless steel inner, anti-acid, anti corrosion characteristics. Corners semicircular is easy to clean, up and down pitch adjustable shelf.</p> <p>5.Electric film heating elements made of carbon fiber ,with advanced cavity preheating, the cavity temperature to reduce temperature.</p> <p>Technical Specifications:</p> <p>Model Number: DNP-9052</p> <p>Volume: 50L</p> <p>Temp. Rang(°C): RT +5 ~ 65°C</p> <p>Temp.Accuracy(°C): ±0.1</p> <p>Temp.Fluctuation(°C): ≤±0.2(in the of range 3~50°C)</p> <p>Timing range: 0 ~ 999 min</p> <p>Voltae (V): 220</p> <p>Power: 0.3KW</p> <p>External Size(mm): 580*530*600</p> <p>Inner Size (mm): 380*380*350</p> <p>Package Type: plywood</p> <p>Single package size: 105X87X82 cm</p> <p>Single gross weight: 70 KG</p> <p>Net. Weight(KG): 31</p> <p>Make: China</p>	



TECHNOCRATIC SERVICES

ISO 9001:2015 Certified Company





TECHNOCRATIC SERVICES

ISO 9001:2015 Certified Company





Advantage silicone seal ring



Thermal motor



Reinforced hinge

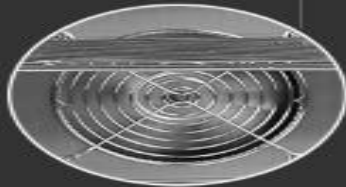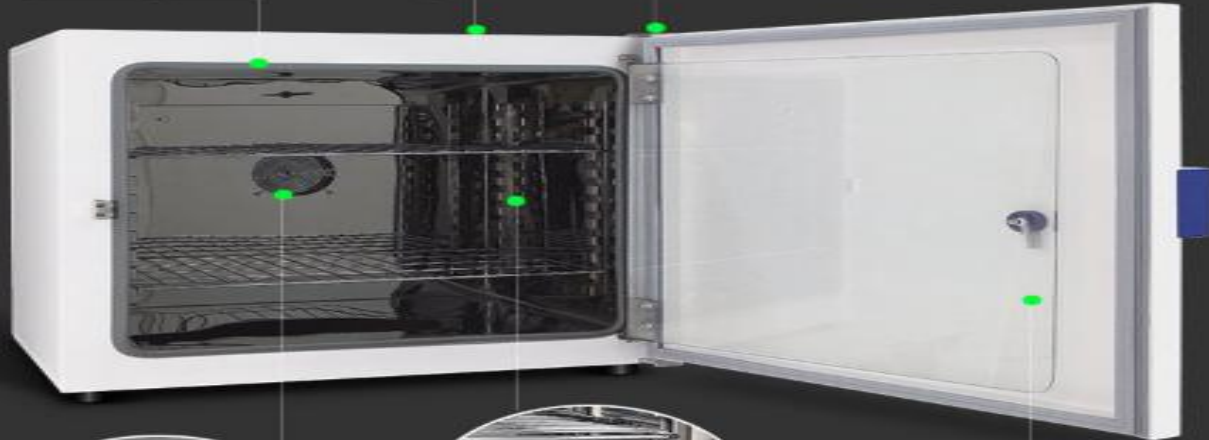


Blower



High quality mirror stainless steel, Spacer spacing adjustable



Tempered glass door, Sealing insulation effect

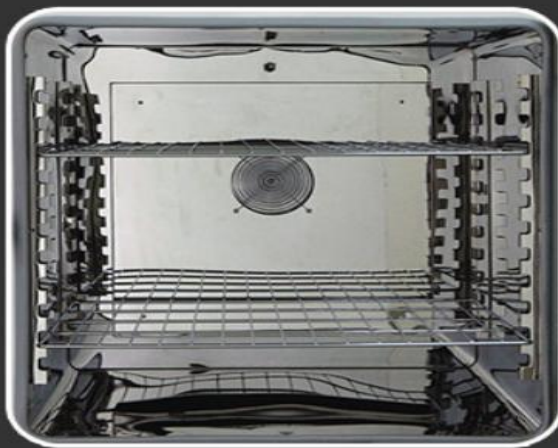

Stainless steel liner

High temperature resistance, fire resistance, corrosion resistance, can be recycled




Double door

Inner tempered glass door, the outer layer of high quality silicone seals, suitable for high temperature culture environment.



TECHNOCRATIC SERVICES

ISO 9001:2015 Certified Company

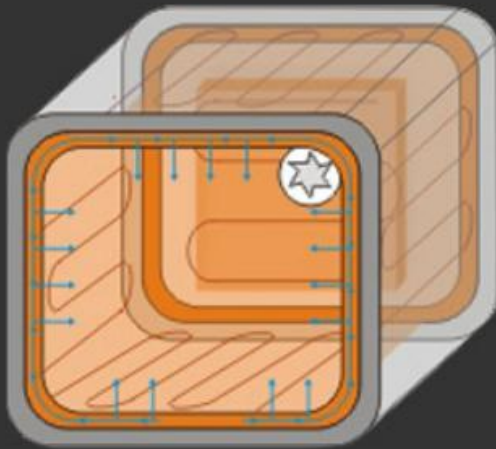

Hot air circulation system

The heating element adopts carbon fiber electric heating film and cooperates with advanced cavity preheating technology to make the inner cavity temperature uniform, reduce temperature fluctuations, and make the temperature of the working chamber uniform.



power cable

Microcomputer interface



TECHNOCRATIC SERVICES

ISO 9001:2015 Certified Company



- ➡ Digital temperature controller (model "C")
- ➡ Intelligent LCD temperature controller

Head Office Karachi

- ☎ +92-21-35001798
- 🌐 www.technocratics.com.pk
- ☎ +92-333-0297007
- ✉ technoratics@gmail.com
- ☎ +92-321-2202241
- ✉ info@technocratics.com.pk
- ☎ +92-300-3972448
- 📍 Plot No. 763, Future Colony Mor, Main Korangi Industrial Area Road, Karachi-75120

Branch Office Lahore

- ☎ 0333-0297007, 0335-4944616, 0321-2202241
- 📍 Plot # 42, Street # 2, Block-C, Asim Town, Harbanspura Underpass Near Ring Road, Lahore-Pakistan.